


Diode laser hair removal — settings & technique

డయోడ్ లేజర్ హేయిర్ రిమూవల్ — సెట్టింగ్స్ & టెక్నిక్

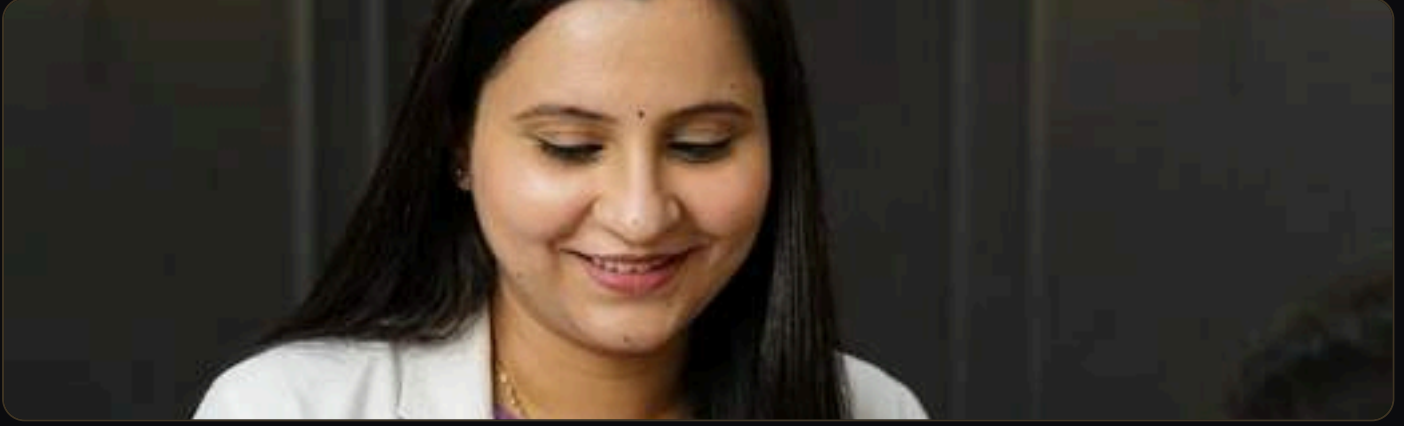

TODAY'S OBJECTIVES · ఈ రోజు లక్ష్యాలు

- ✦ Select safe starting parameters by Fitzpatrick type and body area
- ✦ Perform correct handpiece technique with cooling
- ✦ Counsel the client on the full session course

THEORY · థియరీ

How it works

808 nm diode light is absorbed by melanin in the hair shaft and follicle, heats it and damages the stem cells in the bulge. Only hair in the anagen (growing) phase responds — that is why 6–8 sessions, 4–6 weeks apart on the face and 6–8 weeks on the body, are needed. Grey, white and very fine vellus hair have no melanin and do not respond.

Choosing parameters

In-motion mode: lower fluence, many passes, comfortable and safer for Fitzpatrick IV–V. Static mode: higher fluence, single pass, faster on coarse hair. Typical starting ranges — face 8–12 J/cm², underarms/bikini 10–14, arms/legs 12–16, always adjusted after a test shot. Higher Fitzpatrick → lower fluence, longer pulse, more cooling. Record every setting.

Preparation and technique

Shave the area (never wax/thread for 4 weeks), clean, dry, mark the grid. Contact cooling at 0–5 °C, handpiece flat and full contact, overlap 10–15%, systematic grid so nothing is missed or double-treated. Ask about pain every pass. Endpoint: perifollicular erythema and mild oedema — that is success, not damage.



Hands-on today

- ✦ Test shot + full underarm treatment on a live client under supervision
- ✦ Practise grid technique and overlap on a marked area
- ✦ Record a complete treatment log with settings



Safety & contraindications

- ✦ No LHR over tattoos, moles, active infection, recent tan, or on isotretinoin (6 months)
- ✦ Stop immediately on grey-white change, blistering or severe pain
- ✦ Goggles for everyone, every shot

KEY POINTS TO REMEMBER · గుర్తుంచుకోవలసినవి

- ✦ 808 nm diode, 6–8 sessions
- ✦ Shave, never wax, before session
- ✦ Endpoint = perifollicular erythema
- ✦ Fitz V: low fluence, long pulse, cool hard



Homework

Make a parameter card for 6 body areas at Fitzpatrick IV and V.



Grooming & professional tip

Count your passes out loud. Systematic operators never leave a missed patch — and missed patches are the #1 LHR complaint.

SELF-CHECK · సెల్ఫ్ చెక్

- ✦ Did I follow hygiene and consent for every client today?
- ✦ Which step did I need help with — and what will I drill tonight?
- ✦ Did I write my records and settings before leaving?

Trainer sign: _____ Trainee sign: _____